

MEDICLUS Dental Solution

designed for dentists

www.mediclus.co.kr



MEDICLUS CO., Ltd., was founded in 2011, as a manufacturer of dental chemical materials based on excellent R&D technology. All of MEDICLUS products have acquired patents for domestic and international certifications, such as CE ISO 13485/9001 and KGMP. We have been managed under the strict quality control system by KFDA and SGS. We are exporting to 40 countries. We will continue to strive forward to develop a wide spectrum of dental materials and excellent products for our customers.

MEDICLUS Co., Ltd



■ Temporary Materials

Once-Fil™ Flow
Ezi-Crown™
Any-Temp™
Any-Temp™ Clear

■ Restorative Materials

Any-Com™
Any-Com™ Universal
Any-Com™ Flow
Any-Com™ Flow Rainbow
Bulk-Flo™
MD-Seal™
Any-Etch™ (HV)
Dentto-Etch™ (LV)
Any-Bond™

■ Endodontic Materials

Any-Seal™ RC
One-Fil™ MTA
Any-Pex™
Any-Paste™
Endo@PreP™ EDTA
Any-Core™
Canal-Seek™

■ Impression Materials

Any-Flex™
(Heavy, Light, Bite, Clear)

■ Whitening & ETC

W-Gel™
Any-Flu™
Dentto-Dam™
THE-SEP™
The-Block Out Resin™ (LC)

■ ETC

Ezi-Light™
Silicone Floss™
Resin Display
Tooth Brush

Once-Fil™ Flow

Light-Cured, temporary filling material

Once-Fil Flow has a great flowability and is easily applied to patients.
The elastic force and strength is great for sealing the margin.
Removable in one piece and provide strong masticatory forces.

| SIMPLE |

Easy to handle and to remove in one piece

| SPEEDY |

Time-saving! No instruments needed for occlusion.

| SAFE |

Excellent sealing ability

Contents

#101 Once-Fil Flow 1.2ml syringe*5ea / Disposable tips



#102 Once-Fil 3g syringe*3ea

#103 Ezi-Fil 1.2ml syringe *5ea / Disposable tips


#Color : Blue & Yellow (Natural)



Application 1 For implant case (Curing time 20 ~ 40 sec.) with Ezi-Fil

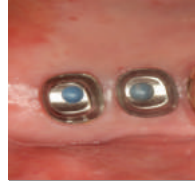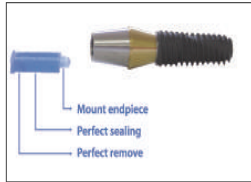

Application 2 For bracket and edge of wire (Curing time 10 ~ 15 sec.)

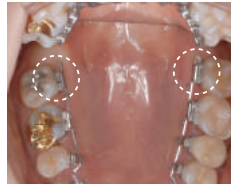

Application 3 For cavity (Curing time 10 ~ 15 sec.)



Ezi-Crown™

Light-Cured, temporary crown and bridge resin

Ezi-Crown LC offers a simple way of making temporary inlay, onlay, crowns and bridges.

Advantages

NO MIXING


Easy and fast

NO STICKY

Excellent handling

NO MESSY

Perfect occlusion accuracy

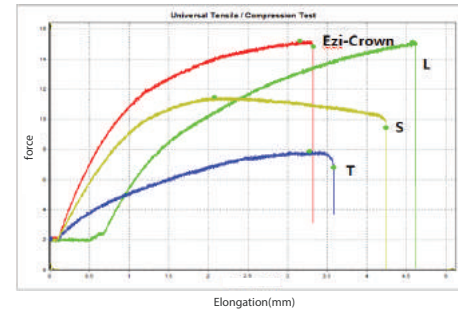

Result



Ezi-Crown 15g*3ea



THE-SEP (Vaseline)




Any-Flex Clear

Application [Single Crown]



1. Apply vaseline on the teeth.
2. Cut off approximate amount of Ezi-Crown.
3. Apply Ezi-Crown.
Press into place, shape and sculpt.

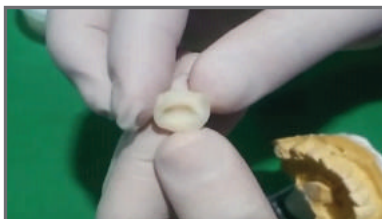

6. Remove the cured Ezi-Crown from the teeth, and apply the curing light again on the inside. (20 sec.)



4. Make the occlusion. Have a patient bite softly to register occlusal surface and adjust



7. If you need polishing and trimming , do that with lab, carbide bur or silicone point.



5. Apply the curing light
(Curing time : 20 sec.)

Contents

#701-1 Kit

Ezi-Crown LC 15g * 3ea

Any-Flex Clear 50ml Cartridge * 1ea

THE-SEP 5ml syringe * 1ea

#701-2 Refill

Ezi-Crown 15g * 3ea

Temporary Materials | 07

Application [Bridge and Crown]



1. Apply Any-Flex Clear on teeth directly or on the tray
(Setting time : 1'30")



2. Prior to preparation or extraction of a teeth, take impression using transparent silicone. Preparation and dry the teeth.


MDCLUS
MEDICLUS Co., Ltd



3. Apply Vaseline (water-soluble separator THE-SEP) on teeth.



4. Inject Ezi-Crown Flow to the prepared impression tray.

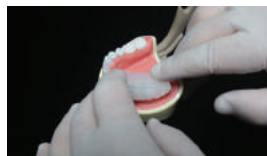


5. Put the impression tray with Ezi-Crown Flow on teeth.



6. Apply curing light
(Curing time : 20 sec.)



7. Removed the cured Ezi-Crown from the teeth and apply the curing light inside again.
(20 sec.)

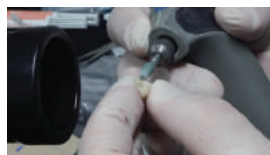

8. If you need polishing and trimming, do with carbide bur or silicone point.



9. Apply temporary cement
(Any-Temp Clear)



THE-SEP (Vaseline)



Any-Temp™

Non-Eugenol, Temporary cement for temporary crowns & bridges.

Advantages

- Easy mix and loss save tip type.
- Enough flowability and easy removal
- Low film thickness
- Setting time $\leq 5\text{min.}$
- Working time $\leq 3\text{min.}$

Contents

- #201-1 Any - Temp NE 10g syringe *2ea , Eco-Tips
- #201-2 Any - Temp NE 10g syringe*3ea , Eco-Tips
- #202 Any - Temp Clear 5ml syringe * 1ea , Auto Mixing tips




Any-TempTM Clear

Temporary cement for implant retained crowns.

Advantages

- Long-term cementation
- Remove it when you need
- Visible in X-ray, good radiopacity
- Easy handling
- Anti-bacterial effect

Application

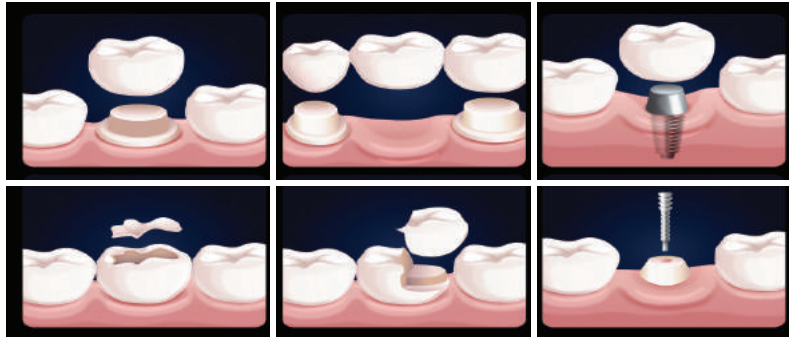

Temporary
Materials

Restorative
Materials

Endodontic
Materials

Impression
Materials

Whitening

ETC



MEDICLUS COMPOSITE TECHNOLOGIES

MEDICLUS composite technique based on a converging technology of **nano-hybrid spherical ceramic filler** and a modified resin monomer, is an excellent resin with outstanding properties in compressive strength, workability, aesthetics and abrasion resistance.



Traditional amorphous Filler

MEDICLUS Filler Particle

Any-Com™

Light-Cured, radiopaque nano-hybrid zirconium composite for the restorative treatment of anterior and posterior teeth.

Application Ranges

- Class I ~ V fillings
- Reconstruction of traumatically affected anterior
- Correction of shape and shade for better anesthetic appearance
- Locking , splinting of loose anterior
- Repairing veneer
- Filling deciduous teeth
- Core build-up under crowns
- Composite inlays
- Depth of cure : 2.5mm (10 seconds)

Contents



#301-1 Kit

- Any-Com 4g syringe * 5ea
- [A1, A2, A3, A3.5, B2 : Shade selectable]
- Any-Etch 3ml syringe * 1ea
- Any-Bond 5ml bottle * 1ea

#301-2 Refill

- Any-Com 4g syringe * 1ea


Physical properties

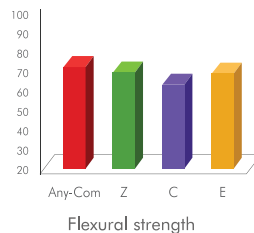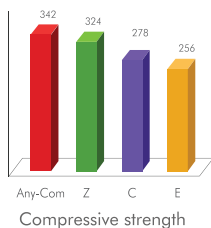

Any-Com™ Universal

Light-Cured, radiopaque nano-hybrid zirconium composite for the restorative treatment of anterior and posterior teeth.

Clinical Reference

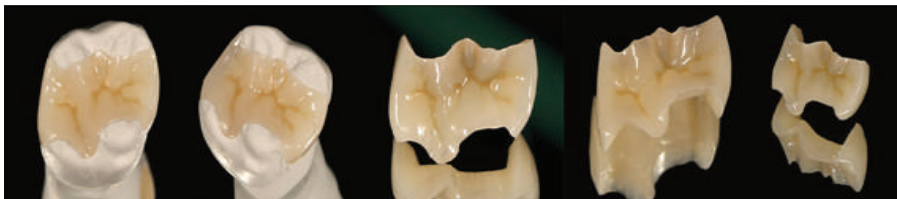

Shade



A1, A2, A3, A3.5, A4, A2E, A3E

B1, B2, B3, B2E, C2, C3

D2, D3, OA2, OA3, TL, WT, OP

Contents

#301-3 Any-Com Universal Kit
 [A1, A2, A3, B2, E1, E2 : Shade selectable]
 Any-Etch 3ml syringe*1ea
 Any-Bond 5ml Bottle * 1ea

Any-Com™ Flow

Light-Cured, radiopaque nano - hybrid flowable zirconium composite for the restorative treatment

Application Ranges

- Restoration for III ~V case
- Reconstruction of traumatically affected anterior
- Base / liner in the Class I , Class II restorations
- Repair of resin, porcelain, and acrylic temporary materials
- Pit and fissure sealant
- Depth of cure : 3mm(10seconds)

Contents

#306-1 Kit

- Any-Com Flow 2g syringe * 4ea
[A1,A2,A3,A3.5 : Shade selectable]
- Any-Etch 1.2ml syringe * 1ea

#306-2 Refill

- Any-Com Flow 2g syringe * 2ea




#306-3 Low Flow

- Any-Com Flow 2g syringe * 2ea
[A1, A2, A3, A3.5, B1, B2]

#306-4 Any-Com Flow Rainbow

- 1g syringe * 5ea
[Pink, Orange, Yellow, Green, Blue, Purple : Shade selectable]

Physical properties

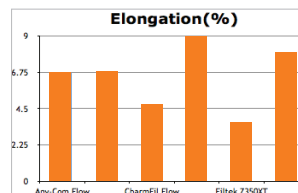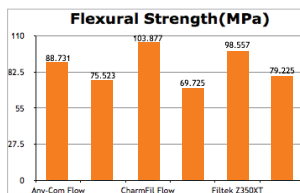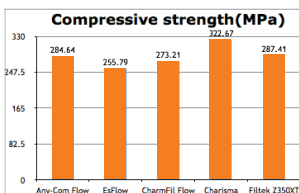

Clinical reference [1]



Clinical reference [2]



Shade



A1, A2, A3, A3.5, A4,
B1, B2, B3, C2, C3
D2, OA2, OA3, TL, OWT, YTL, RTL



Rainbow Shade




Bulk-Flo™

Bulk fill flowable resin



It can be used as a base-liner for both direct and indirect restoration up to 4mm a single layer placement.
It is also suitable as a liner for deep cavity.(Class cavity)

- ✓ Resin-based
- ✓ Light-cure
- ✓ For bases and liners
- ✓ For both anterior and posterior use
- ✓ Radiopaque
- ✓ Flowable


MEDICLUS technology

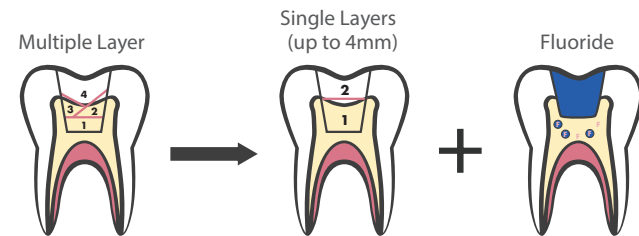

Results [Core Application]

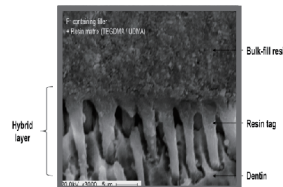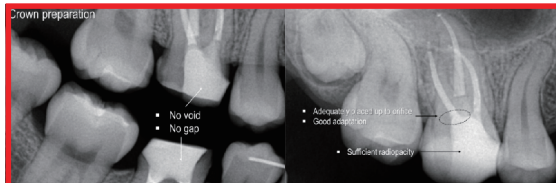


Showing resin tag

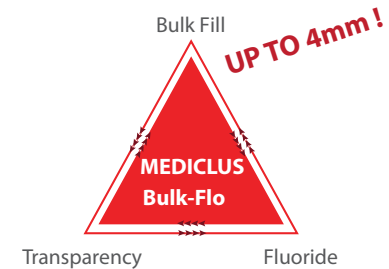


Applications

Step1

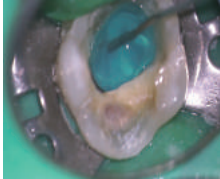


Apply Etching Gel

Step2

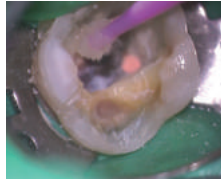

Apply Primer and Bond

Step3

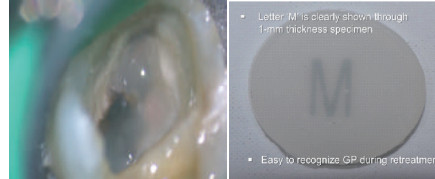

Transparency shade

Step4



Apply Bulk-Flo
(Universal)

■ For Direct restoration with a composite

1. Place any choice of composite restorative material over the Bulk-Flo any light-cure according to the manufactures 's instruction.
2. At least 2mm of the composite surface restorative material is required.

■ For Indirect restoration with a prosthesis

1. Prepare the tooth in the usual manner.
2. Take an impression with agar; alginate combination or silicone-base impression material to fabricate a prosthesis.
3. Temporize with a temporary material until cementation of the prosthesis the manufactures 's instructions.
(Avoid using eugenol-containing temporary material)

Contents

#307 Bulk-Flo 2g syringe * 3ea.
Disposable tips



[Shade]

A1, A2 , A3 , Blue , Universal


MD-Seal™

Light cured, pit& fissure sealant


Advantage

- ✓ Fast and safe application
- ✓ Low polymerization shrinkage
- ✓ High bonding strength
- ✓ Appropriate viscosity
- ✓ Nano filler
- ✓ Opaque shade

Physical properties

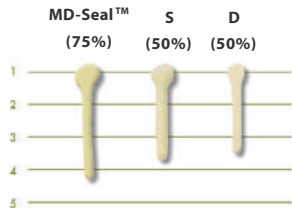

Resin content and flowability

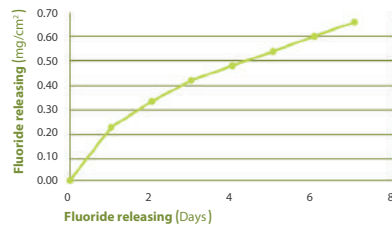

Opacity



Contents

#308 MD-Seal 2g syringe * 2ea
Disposable tips



Results

<Result1>



<Result2>





Any-Etch™ (HV)

37% Phosphoric acid etching agent with unique optimum viscosity

Advantages

- A cream-typed etchant preventing corrosion
- Forming the best etching depth(Dentin 10sec./ Enamel 15sec.)
- Available for both in dentin and enamel
- Allows precise placement in any location

Contents

- #304-1 Any-Etch(HV)1.2ml syringe * 5ea
- #304-2 Any-Etch(HV)3ml syringe * 5ea
- #305 Dentto-Etch 3ml syringe * 5ea



Dentto-Etch



Any-Etch(HV)
[SELECT- ETCH]



Select Etching



Allows precise placement
in any location

"SELECT ETCH"



Flow Down ?? NO!!



Perfect Viscosity!!!

Any-Bond™

5th generation bonding system one-step bonding agent for priming and bonding

Advantages

- Bonding for direct composite restoration
- Bonding for porcelain, composite, metal repair, setting amalgam and post
- Bonding for porcelain veneer

Contents

#303 Any-Bond 5ml Bottle * 1ea



Temporary
Materials

Restorative
Materials

Endodontic
Materials

Impression
Materials

Whitening





ETC

Any-Seal™ RC

Resin based root canal sealer

Advantages




- Excellent sealing ability
- Excellent biocompatibility
- Non-staining to the resin
- Good radiopacity
- Easy and fast handling
- Working time 1hr ~ 1hr 30min
- Setting time : 3hr ~ 4hr at 37°C

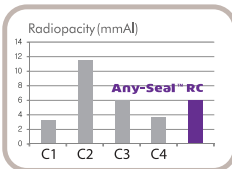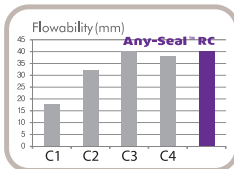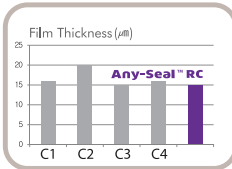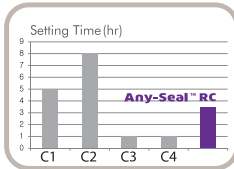

Contents

#503 Any-Seal RC 10g syringe * 1ea
Disposable tips

Application

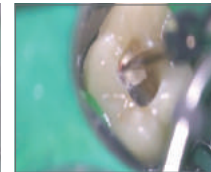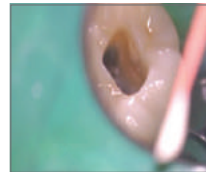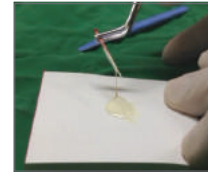

One-Fil™ MTA

Pre-mixed root canal filling material

One-Fil MTA is a paste-type root canal sealer based on Calcium silicate cement that has excellent physical and biological properties

Advantages

- Simple premixed injectable paste type
- White color
- Setting time : within 2 hours
- High radiopacity : 9mm / Al
- Adequate flowability : 21mm(by ISO6876)
- Antibacterial effect
- Dimensional stability
- Biocompatibility
- Outstanding sealing property
- Promoting hard-tissue formation
- Insolubility

Contents

#508 One-Fil MTA 2.0g syringe * 1ea
Disposable tips

Temporary
MaterialsRestorative
MaterialsEndodontic
MaterialsImpression
Materials

Whitening

ETC

Any-Pex™

Oil-Base Calcium hydroxide paste with Iodoform temporary root canal filling materials.

Advantages

- Excellent antibacterial effect
- Excellent radiopacity
- Premixed paste and syringe type

Applications

- Root canal disinfection
- Apexification
- Detection of canal aberrations

Contents

#504 Any-Pex 2.2g syringe * 1ea
Disposable tips



Any-Paste™

Water-based Calcium hydroxide paste with Barium Sulfate temporary root canal filling materials

Advantages

- Excellent antibacterial effect
- Excellent radiopacity
- Premixed paste and syringe type

Applications

- Root canal disinfection
- Apexification
- Detection of canal aberrations

Contents

#505 Any-Paste 2g syringe * 1 ea
Disposable tips



Temporary
Materials

Restorative
Materials

Endodontic
Materials

Impression
Materials

Whitening

ETC

Endo@PreP™ EDTA

EDTA-based root canal cleanser

For root canal cleaning and preparation


A patented root canal cleaning and preparation gel with a great chelation and lubricating effect.

Dental cleanser for removal of dentulpulp tissues in root canal treatment.


Advantages

- Removal of inorganic substances by chelation(EDTA-Ca)
- Easy application with lubricating action
- Provide a clean, smear-free canal wall


Application

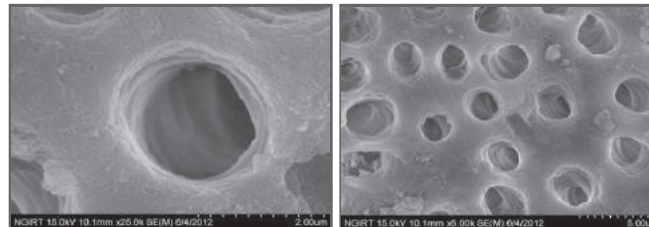

①Apply Endo@prep Gel

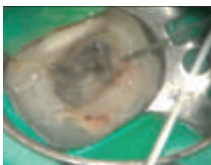

②Apply NaOCl

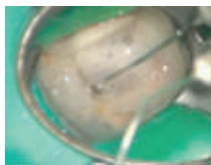

③Bubbling Effect



[Reference]

Kim HJ, Park SJ, Park SH, Hwang YC, Yu MK, Min KS.
Efficacy of flowable gel-type EDTA at removing the smear layer
and inorganic debris under manual dynamic activation.
J Endod 2013;39:910-4. By Prof. Kyung-San Min in Chonbuk
National University school of Dentistry.

Basic Research—Technology

Efficacy of Flowable Gel-type EDTA at Removing the Smear Layer and Inorganic Debris under Manual Dynamic Activation

Hee-Jin Kim, DDS, MSD,* Su-Jung Park, DDS, PhD,* Sang-Hyuk Park, DDS, PhD,^{1,2}
Yun-Chan Hwang, DDS, PhD,^{3,1} Mi-Kyung Yu, DDS, PhD,⁵ and Kyung-San Min, DDS, PhD⁵

Abstract

Introduction: A flowable gel-type EDTA solution containing urea peroxide and polyacrylic acid was recently introduced into the endodontic market. However, its efficacy for removing the smear layer and inorganic debris remains unknown. This study was performed to investigate the relative efficacies of the flowable gel-type and liquid-type EDTA solutions for removal of the smear layer and inorganic debris. We also evaluated the effects of manual dynamic activation (MDA). **Methods:** Wettability was evaluated by measuring the contact angle. The incidence of accidental extrusion of irrigant was determined. The effervescent effect was evaluated by mixing the solutions with sodium hypochlorite. The efficacies of the EDTA solutions at removing the smear layer and inorganic debris were evaluated by scanning electron microscopic examination. **Results:** The contact angles of the 2 EDTA solutions did not differ significantly throughout the experiment ($P > .05$). Accidental extrusion occurred 4 times for the liquid-type EDTA but never for the gel-type EDTA. The gel-type but not the liquid-type EDTA showed an effervescent effect. The EDTA/MDA treatment combinations did not produce significantly different smear layer scores ($P > .05$). However, the debris scores for the coronal and middle parts were significantly lower for the gel-type EDTA with MDA than for the liquid-type EDTA without MDA ($P < .05$). **Conclusions:** Our results suggest that the newly introduced gel-type EDTA might be an acceptable irrigant for removing the smear layer and inorganic debris present on the root canal wall. (J Endod 2013;39:910-914)

From the *Department of Conservative Dentistry, Wonkwang University School of Dentistry, Gwangju, Korea; †Department of Conservative Dentistry, Kyunghee University College of Dentistry, and Yonsei Woon Hae Medical Center, Seoul, Korea; ‡Department of Conservative Dentistry, School of Dentistry, and Dental Science Research Institute, Chonnam National University, Gwangju, Korea; §Department of Conservative Dentistry, Chonbuk National University School of Dentistry, Jeonju, Korea; ||In Jin Kim and Su Jung Park contributed equally as first authors.

Address requests for reprints to Dr Kyung-San Min, Department of Conservative Dentistry, Chonbuk National University School of Dentistry, 567 Baekje-daero, Daejeon-shi, Jeonju-shi, Jeollabuk-do 561-758, Republic of Korea. E-mail address: endodontid@gmail.com
0099-2398/\$ - see front matter
Copyright © 2013 American Association of Endodontists.
http://dx.doi.org/10.1016/j.jen.2013.04.018

Key Words

Debris, EDTA, effervescent, gel, manual dynamic activation, smear layer

Chelating agents were introduced into endodontics to facilitate the preparation of narrow and calcified root canals. Of these products, a liquid form of EDTA is the most widely used as a final irrigant for removing the smear layer and other inorganic substances (1). Although there is some controversy regarding the removal of the smear layer overlying dentinal tubules (2), previous studies have shown that the smear layer adversely affects disinfection and increases microleakage after canal obturation (3-5). Many liquid EDTA solutions have been introduced for endodontic use. However, most of them should be transferred to a syringe before use. In addition, liquid-type solutions as a factor of their lack of viscosity may be extruded unexpectedly if the clinician inadvertently depresses the plunger of the syringe. A very low-viscosity, flowable gel-type EDTA (Endo-Prep EDTA Gel; Medialis, Chungju, Korea) containing urea peroxide and polyacrylic acid, a so-called carbonizer, was recently introduced. This gel-type EDTA is packaged in a preloaded syringe and can be used with various sizes of needles. Polyacrylic acid is added as a thickener to provide viscosity and a lubrication effect. Furthermore, the urea peroxide in this formulation is known to produce effervescence if it contacts sodium hypochlorite (NaOCl) (6). The effervescent reaction can push debris out of the root canal through the path of least resistance into the chamber (7). However, more viscous EDTA solutions may fail to reach the apical areas of the root canal when used as irrigants. Therefore, most viscous chelators are used predominantly to aid the negotiation of calcified canals or decrease frictional stress during root canal preparation using rotary instruments, and the viscosity of any solution to be used as an irrigant for smear layer removal should be regulated carefully. However, there has been no study showing the efficacy of the flowable gel-type EDTA at removing the smear layer and inorganic debris.

An irrigant must directly contact the canal walls in order to be effective. However, it is often difficult to ensure that the irrigant reaches the apical portion of the canal because of the so-called vapor lock effect (8, 9). Research has shown that gently moving a well-fitting gutta-percha master cone up and down in short 2- to 3-mm strokes (manual dynamic activation [MDA]) within an instrumented canal can produce an effective hydrodynamic effect and significantly improve the displacement and exchange of any desired reagent (10, 11). Although MDA has been advocated as a method for canal irrigation because of its simplicity and cost-effectiveness, there is little information about its effects on the smear layer and inorganic debris when used with the gel-type EDTA. This study aimed to investigate the efficacy of the newly developed flowable gel-type EDTA versus that of the conventional 17% liquid-type EDTA solution at removing the smear layer and inorganic debris. We further explored whether MDA could improve the efficacy of smear layer and debris removal.

Materials and Methods

Dentin Specimen Fabrication

Freshly extracted human third molars without fractures or carious lesions were used. The teeth were embedded using autopolymerizing acrylic resin (Ortho-Jet; Ling Dental, Wheeling, IL). The specimens were then polished under distilled water



Contents

- #501 Gel Type 9g syringe * 2ea / Disposable tip
- #502 Cream Type 9g syringe * 2ea / Disposable tips

Any-Core™

Any-Core dual syringe
Dual Cured, core build up composite resin

Advantages


- Optimum viscosity and flowability
- Dual-cured, Light-cured and self-cured
- Radiopaque
- Perfect occlusion accuracy
- Working time : 2min
- Light curing time : 20sec.(2mm)/40sec.(4mm)
- Chemical curing time : 5min(Intra oral)
- Compressive strength : 338Mpa
- Flexual strength : 167Mpa
- Various shade : [Blue, Tooth, White]



Contents

- #801 Kit
 - Any-Core 9g syringe * 2ea
 - Disposable tips
- #802 Refill
 - Any-Core 9g syringe * 1ea/No tips

Clinical Reference

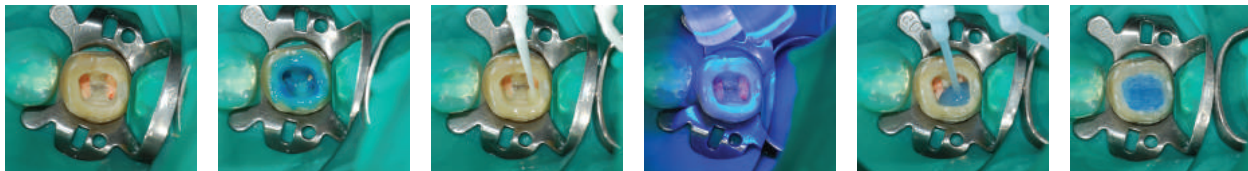

Canal-Seek™

Canal Orifice Indicator

Blue dye contained in the product penetrates the root canal orifices and dyes them blue which makes it easier to find the main root canal and detect the cracks.

How to use

- 1.Remove dentine where the root canal orifices are expected.(picture1)
- 2.Insert 1-3 drops into the bottom of the pulp cavity,wait 5-10 seconds.(Picture2)
- 3.Rinse and dry the pulp cavity(Picture3)
- 4.The areas which show up in blue indicate the root canal orifices.(Picture4)
- 5.After use pull back the pluger of the syringe to prevent overflow of product, remove applicator,close the syringe with cap.

Contents

#507 Canal-Seek 1.2ml syringe * 2ea
Disposable tips



Application



Any-Flex™ (Heavy, Light, Bite, Clear)

Any-Flex are hydrophilic vinyl polysiloxane impression materials in the types of Heavy, Light, Bite and Clear.



Contents

- #901 Heavy-body 50ml Cartridge * 4ea
- #902 Light-body 50ml Cartridge * 4ea
- #903 Bite 50ml Cartridge * 4ea
- #904 Clear 50ml Cartridge * 4ea

Time of application

Items	Working Time (23℃)	Setting Time (Intra-oral, 23℃)	Shore A
Light	1:30	4:30	65
Heavy	1:30	4:30	61
Bite	0:30	0:60	40
Clear	0:30	1:00	70

[Light]

It is used for the one-step(Light body and Heavy body)or two-step(Light body and Putty)impression taking of individual or whole teeth in the mouth.

[Heavy]

It is used for the one-step taking of whole teeth or individual tooth with Light body.
And it is used with Light-body to make a precise impression.

[Bite]

It is a bite registration impression material to measure of the occlusal surface, impression of the teeth of a 3-dimension position for the maxilla and manible.It has a short polymerization time and a high final handress. It is easy to use and removal of the low elastic deformation is suitable for bite refistration.

[Clear]

It is used for one-step taking of serveral teeth.Own characteristic allows penetration of dental curing light's source, so it is possible to cure the inside material.

W-Gel™

W-Gel is designed as multi-function gel

Advantages

- Maximum strength whitening gel(Apply 10min.*3times)
- Polishing agent with rubber cups
- Sticky & high viscous gel
- Orange flavor

Any-Flu™

Fluoride gel with potassium nitrate

Any-Flu can eliminate the sensitivity from tooth brush abrasion, chemical changes and tooth bleaching.

Advantages

- Immediate result
- Eliminate sensitivity
- Sticky and high viscous gel
- Mint flavor



Contents

#401 W-Gel 5ml syringe * 2ea
Mixing bowl
Mixing brush

Potassium Nitrate
3%



Contents

#403 Any-Flu 1.2ml syringe * 5ea
Disposable Tips

Dentto-Dam™

Light-cured, Liquid Rubber dam

Advantages

- Dentto-Dam is a light cured resin which protects gums for in-office bleaching and endodontic treatment
- Gingival protection for whitening, etching micro-abrasion and sandblast etc.
- Perfect sealing for rubber dam
- Endo - Wall

Applications



Contents

- #402 Dentto-Dam 1.2ml syringe*5ea
- Disposable tips



Temporary
Materials

Restorative
Materials

Endodontic
Materials

Impression
Materials

Whitening

ETC


THE-SEP™

Separating agent and oxygen-blockers working together

Advantage

- Water-soluble Vaseline/Separator
- Easy removal of excess cement
- A temporary prosthesis can be used for isolation or oral agent such as laminate
- An implant and the adhesive strength of the temporary adjustment cement
- It can be used as an oxygen-blocker

Applications

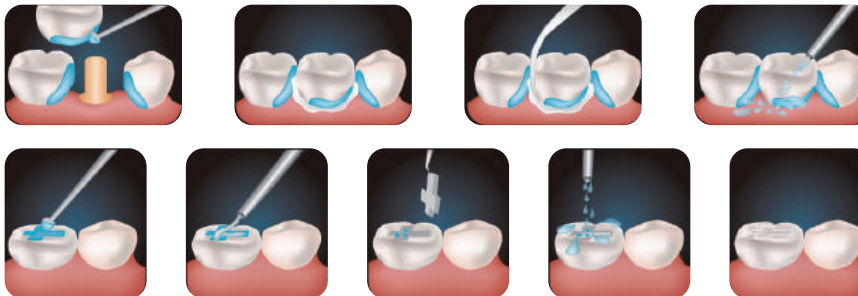

Contents

#701 THE-SEP 5ml syringe*2ea
[shade : Blue , White]



The-Block Out Resin™ (LC)

The-Block out resin is a light cured resin, which was originally developed as a rubber dam substitute. It can be used to make the bleaching tray space, undercut and shaping model for a clear aligner on the model.

Advantages

- Excellent viscosity (Non-runny)
- Provides reservoir space for bleaching tray
- Use for other laboratory procedures
- Blue color

Contents

#601 The-Block Out Resin 1.2ml syringe * 5ea
Disposable Tipe



For Under Cut



For Clear Aligner



For Bleaching Tray

Ezi-Light™

Ezi-Light is a high power translux LED used oral caries detector.

Advantages

- 3W High Power Laser LED
- High concentration LENS
- Wire less & Chargeable Battery
- Easy to find out tooth problems with X-ray irradiation
- Simple Use

Applications

- Caries Detecting
- Root canal pathway exam
- Fractured tooth or weakened root exam
- Caculi and overflow Cement exam
- Implant marginal exam
- All Ceramic& Zirconium fractured exam


Clinical Reference



Root Canal Pathway exam



Under gingival calculus detecting

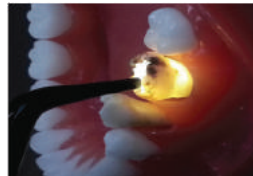

Caries Detecting



Fractured Teeth exam



Contents

#M01 Ezi-Light 1set

Silicon Dental Floss™


Innovative Personal Scaling System SILICONE DENTAL FLOSS

Advantages

| SAFETY | Ergonomically-designed floss is made with 100% silicone material. It does not damage gums, making it easy to use.

| DESIGN | Rounded-shaped handle allows beginners to apply it easily.

| HYGIENE | Made with 100% silicone, it is non-toxic and does not harm the human body. It is also antibacterial and hygienic. Wash it with running water after flossing for more hygienic use.

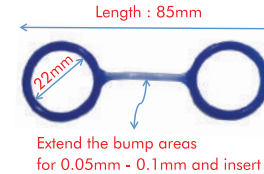

Resin Display



Tooth Brush

Temporary
MaterialsRestorative
MaterialsEndodontic
MaterialsImpression
Materials

Whitening

ETC

No.1210 , 134 Gongdan-ro , Heungdeok-gu , Cheongju-si , Chungbuk , Korea
TEL. +82-43-211-2877 FAX. +82-43-211-2866
E-mail. sales@mdclus.com